



OMNISENSE 7000P

MEETING THE CHALLENGE OF PEDIATRIC BONE HEALTH

 **Sunlight**

Technology solutions at the point of care



OMNISENSE 7000P

BUILDING A HEALTHY FUTURE

- Radiation-free
- Quick and child-friendly testing
- High-precision measurement
- Monitors bone development
- Sensitive to lifestyle-related changes in the bone

BONE HEALTH FOR A LIFETIME


Proper skeletal development during childhood and adolescence is a key determinant of bone health throughout life. Even among otherwise healthy young people, inadequate physical activity and poor nutrition can compromise bone development and lead to fractures¹, poor skeletal health and osteoporosis later in life.

Omnisense 7000P is designed for the management of bone development in children and adolescents. It is the only solution available today for routine children's bone assessment. Omnisense is radiation-free, quick and easy to use, and provides physicians with an essential tool to monitor children's bone development.

REACHING THE PEAK

Reaching maximal peak bone strength by early adulthood minimizes the risk of osteoporosis later in life². Building bone towards this peak starts in childhood, with an active lifestyle and a balanced diet with adequate levels of calcium.³⁻⁵

At-risk children have even lower chances of attaining a maximal peak bone strength. Risk factors include premature birth, a low birth weight, obesity, delayed sexual development, and a variety of medical disorders.⁶⁻⁷

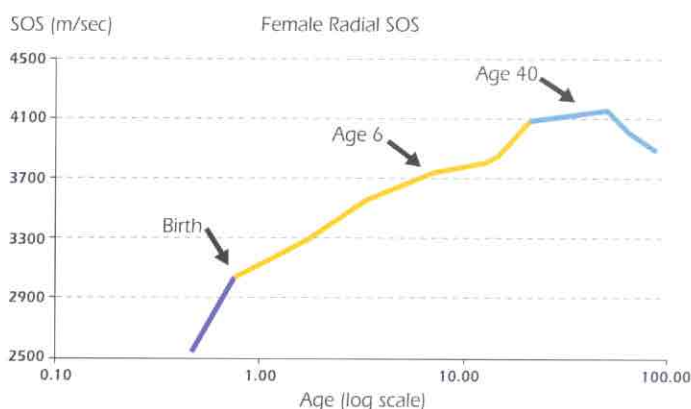

The normal bone growth curve shows a sharp increase in bone strength throughout childhood, culminating in peak bone strength in the late 30s and decreasing thereafter.


MONITORING BONE STRENGTH

Sunlight Omnisense™ 7000P was developed to meet the increasing need for routine pediatric bone health assessment. Omnisense measurement is sensitive to different levels of calcium intake and physical activity. This assessment can help physicians track the development of bone and decide on diet and lifestyle changes when necessary. The precision of Omnisense enables regular monitoring of pediatric bone strength development.


MEETING THE NEED

Omnisense measurement uses Omnipath™, a patented axial transmission technology that measures the speed of sound of ultrasound waves propagating along the bone. Omnipath™ technology eliminates the effects of soft tissue, enabling accurate measurement. This measurement is independent of body size and bone size.

Omnisense boasts a unique pediatric database that compares each result with age- and gender-matched average values, and quickly generates an easy-to-read report.

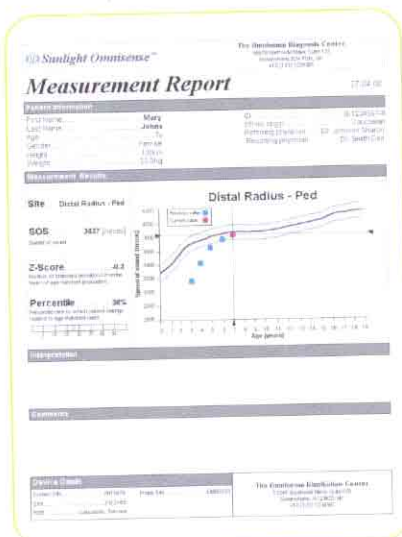

A detailed measurement report portrays results in SOS, Z-score, and percentile formats and shows comparisons with previous measurements.

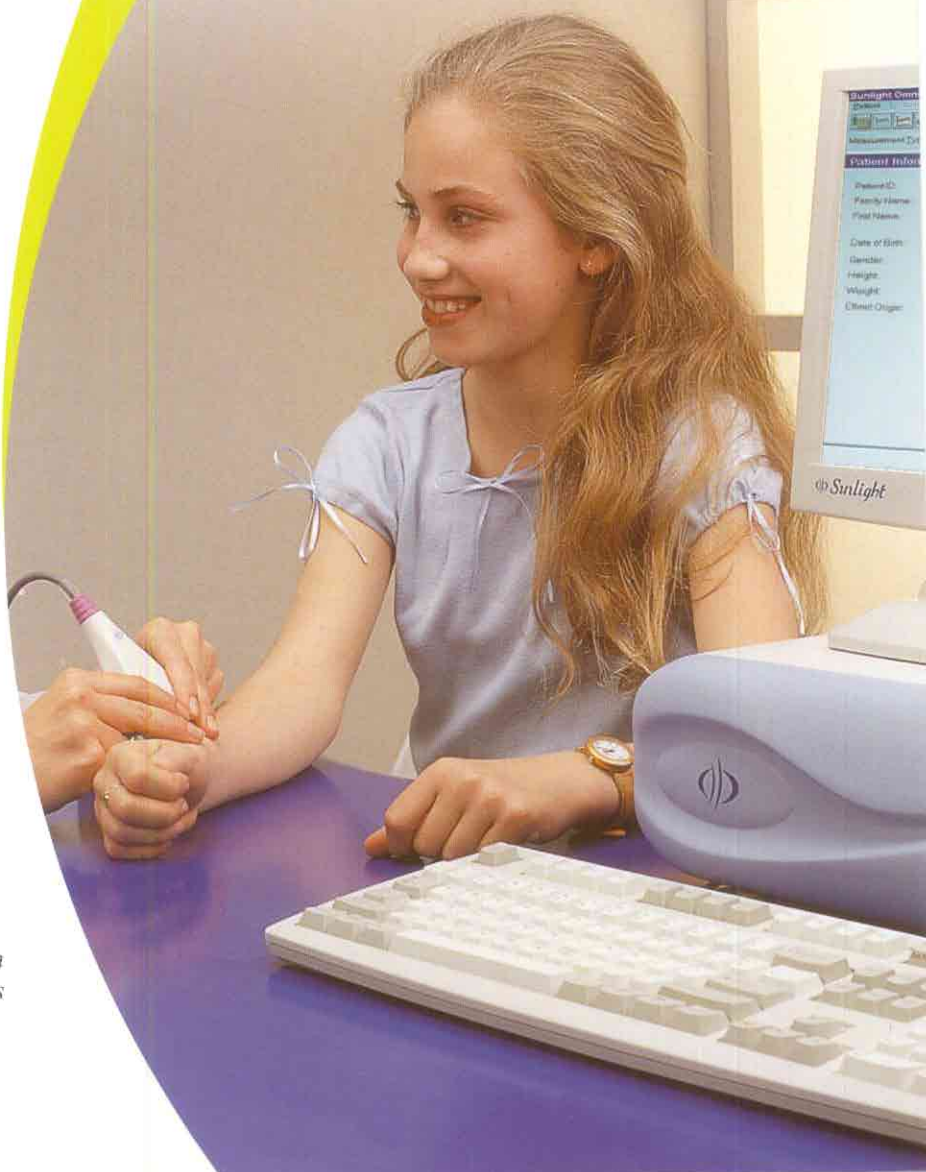

SPECIFICATIONS

Model	Sunlight Omnisense™ 7000P
Precision	RMS CV = 0.32% - 0.36% in-vivo precision, depending on site
Measurement sites	Distal 1/3 radius (forearm); Mid-shaft tibia (lower leg)
Technology	Quantitative ultrasound, using Omnipath™ axial transmission technology
Measured parameter	Speed of sound (SOS), expressed in m/sec
Scan time	Less than one minute per skeletal site
Data analysis	Compares SOS results with reference database and reports Z-scores and percentiles.
Display	Flat color display monitor (LCD TFT)
Power	115 or 230V autoswitchable; ~50-60Hz
Power consumption	Approximately 85 VA (including LCD 80 VA (system alone)
Unit dimensions	39cm x 13cm x 33cm
Main unit weight	Approximately 7 kg

Specifications subject to change without prior notice. Sunlight Omnisense is a trademark of Sunlight Medical Ltd. Other product names are trademarks of their respective holders.

ABOUT SUNLIGHT

Sunlight Medical Ltd. is an international company with offices in the USA, Germany, Israel, and China. Sunlight develops, manufactures and markets advanced, high technology medical devices for the use in primary care facilities and hospitals.

All Sunlight products are tested for safety and EMC requirements applicable to medical devices and comply with both international (IEC 601, CE MDD) and FDA requirements for medical devices.

REFERENCES

- 1 Bailey D.A., J.H. Wedge, R.G. McCulloch, A.D. Martin, S.C. Bernhardtson, Epidemiology of Fractures of the Distal End of the Radius in Children as Associated with Growth, *Journal of Bone and Joint Surgery (Am.)*, 1989, 71, 1225-312.
- 2 NIH Consensus Development Conference on Osteoporosis Prevention, Diagnosis, and Therapy, March 27-29, 2000, 2.
- 3 Calcium Requirements of Infants, Children, and Adolescents (RE9904), *American Academy of Pediatrics, Pediatrics*, Volume 104, No. 5, November 1999, pp. 1152-1157.
- 4 Physical Activity and Health - A Report of the U.S. Surgeon General, 1996.
- 5 Mark S., H. Link, Reducing Osteoporosis: Prevention During Childhood and Adolescence, *Bulletin of World Health Organization*, 1999, 77:423-42.
- 6 Steelman L., P. Zeitler, Osteoporosis in Pediatrics, *Pediatrics in Review*, 2001, 22, 56-65, 2.
- 7 Elakim A., D. Nemet, B. Wolach, Quantitative Ultrasound Measurements of Bone Strength in Obese Children and Adolescents, *Journal of Pediatric Endocrinology & Metabolism*, 2001 Feb; 14(2): 159-642.



Sunlight Medical Ltd.
5 Tuval Street
Tel Aviv, 61251, Israel
Tel: (972)3-5620606
Fax: (972)3-5620607
info@sunlightnet.com

Sunlight Medical Inc.
100 Davidson Avenue, Suite 108
Somerset, NJ 08873, USA
Tel: (1)732-560-8770
Fax: (1)732-560-9462
Toll Free: 1-800-750-6011
info@sunlightnet.com

Sunlight Medical (Deutschland) GmbH
Grüner-Turm-Strasse 16, D-88212
Ravensburg, Germany
Fon: (49)751-363-6230
Fax: (49)751-363-6239
Toll Free: 0800-974-9099
info.de@sunlightnet.com

Sunlight Medical Asia
Shanghai Representative Office
Room 1984, Building B, 100 ZunYi Road
Shanghai 200051, China
Tel: 86-21-62370508
Fax: 86-21-62370509
infoasia@sunlightnet.com

w w w . s u n l i g h t m e d i c a l . c o m

Product is available for research use.
Not available for sale in the U.S pending FDA registration.